



SENTIVA™

This changes everything.

Treating patients with drug-resistant
epilepsy just got easier.



**Smart Technology.
Customisable Therapy.
Proven Results.¹**

SENTIVA™


SenTiva is the smallest and lightest responsive therapy for treatment of drug-resistant epilepsy



Smart detection. Rapid response.

More seizure-free moments.

SenTiva with AutoStim responds to heart rate increases that may be associated with seizures²

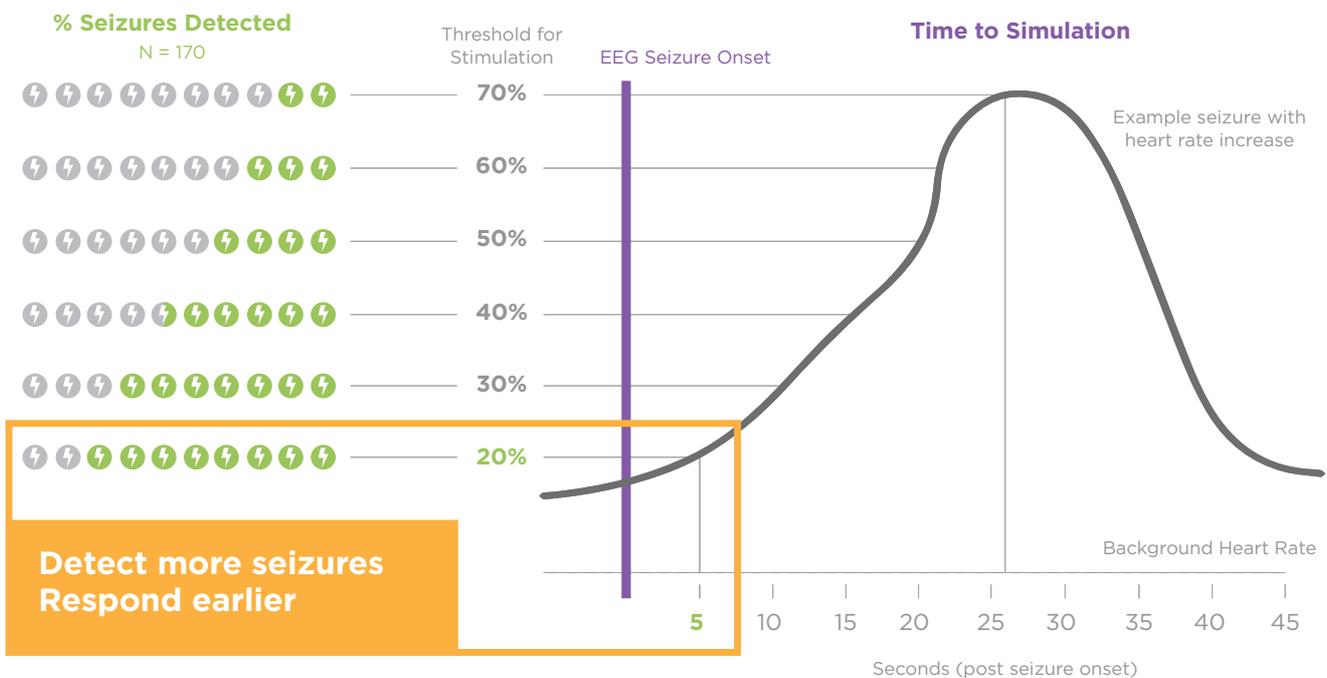

At a **20% threshold**, AutoStim detected **80% of seizures** and stimulated **5 seconds post seizure onset** (median latency)³

Potential to stop or shorten a seizure

Over 60% of treated seizures (N=28/46) ended during automatic stimulation³


Stimulation-associated desynchronisation of focal seizure*³

Case 1⁴

Age: 27

History:

- Bitemporal symptomatic epilepsy with frequent focal seizures without aura

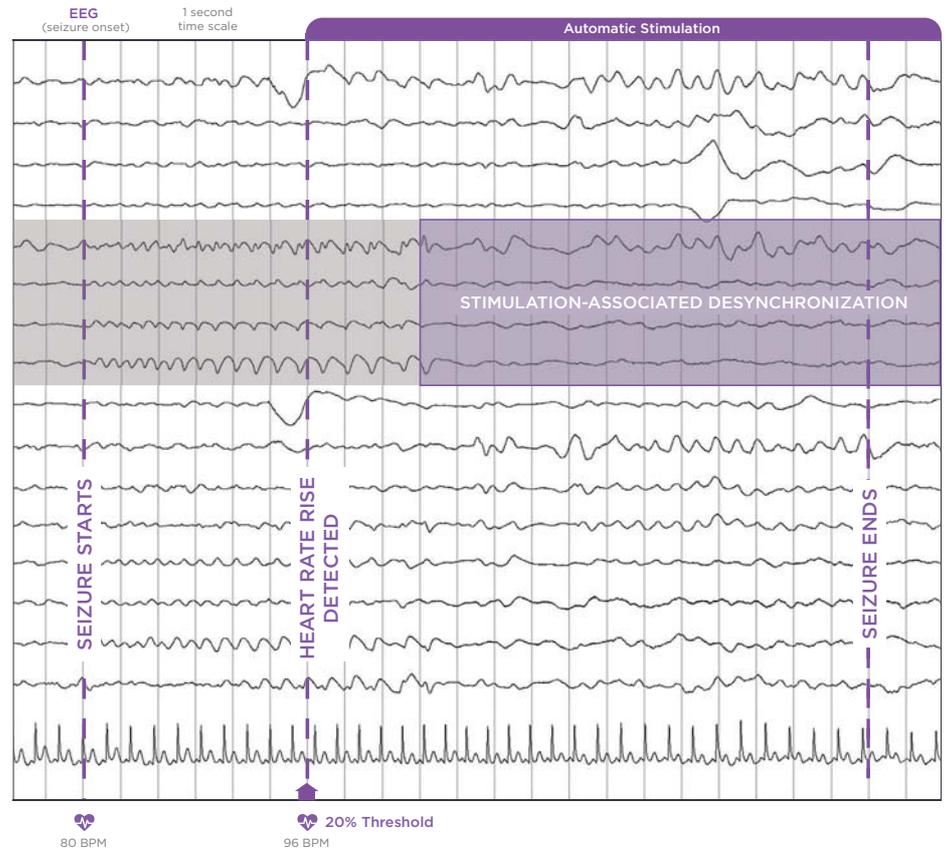

Threshold choice affects timing of stimulation*⁴

Earlier stimulation correlates with shorter seizures³


Case 2⁴

Age: 24

Age at Epilepsy Onset: 16


History:

- Frequent Focal to Tonic-Clonic Seizures
- Bihemispheric Epilepsy of Unknown Origin

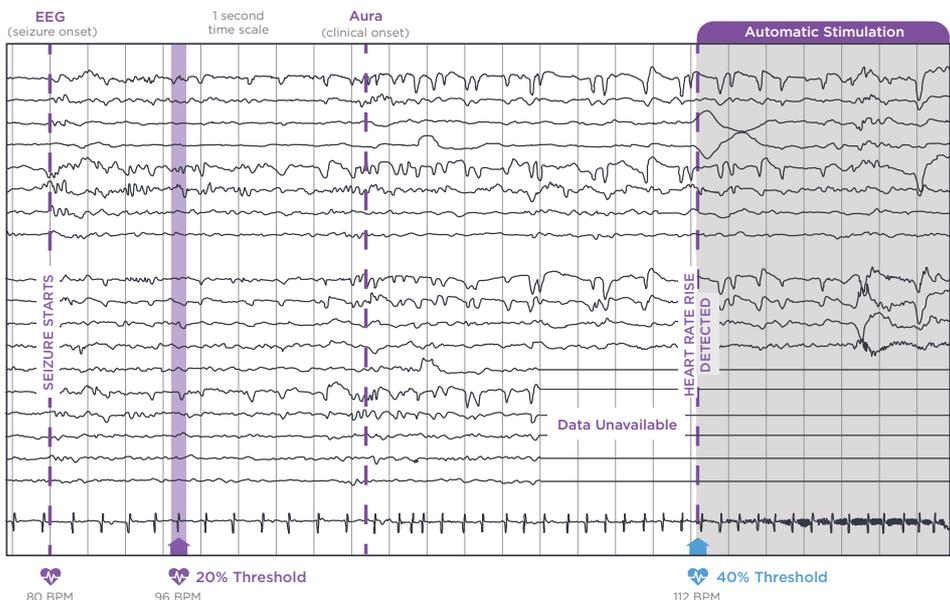

*VNS Therapy with AutoStim
 **Individual study patient example. Results may vary.

Confidently deliver VNS Therapy with one touch²

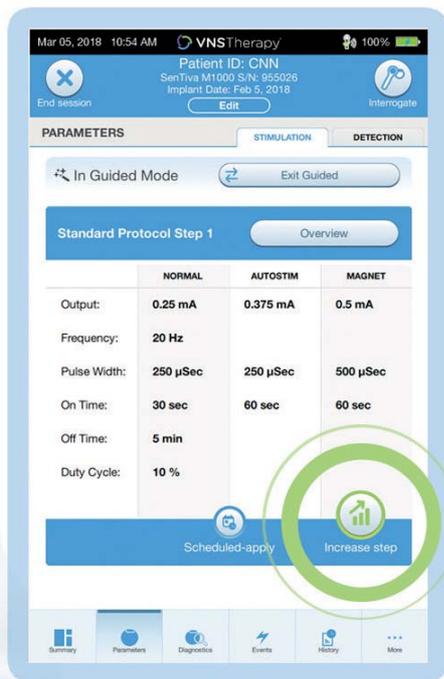

Guided Programming

- Based on an CE-marked approved protocol
- Safe and easy pathway towards achieving targeted therapy levels through one-touch programming
 - ▶ Parameter defaults chosen to minimise side effects



Treating patients with drug-resistant epilepsy just got easier


Safely titrate multiple steps without office visits²

Scheduled Programming

- › Optional feature allowing the healthcare provider to program a titration schedule in one office visit
- › Titration changes will be applied at the defined intervals while the patient lives their life

? Do your patients:


- › Miss dosing appointments?
- › Have difficulties with travel to appointments?

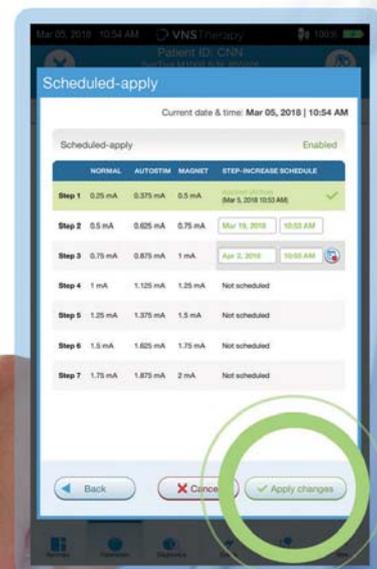

Quickly visualise patient trends²

Events & Trends

- › Quickly view data at follow-up visits on events that may be associated with seizures
- › Easily export patient data and trends

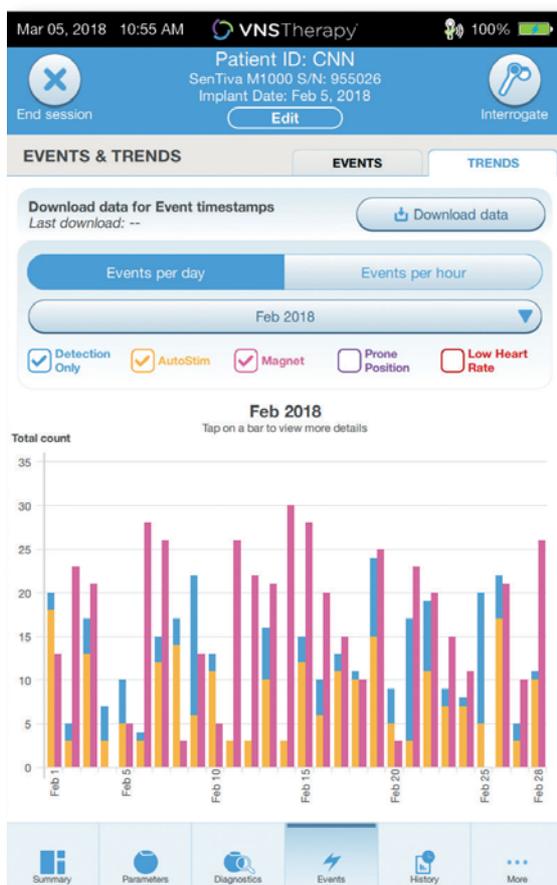

Generator shown is actual size


Customised control for when your patients need it²

Day and Night Programming

- › Two independent sets of parameters can be customised and will alternate based on time of day
- › Select which parameters will change, what time the change will occur, and the duration of the change



› Which patients could benefit from this customised control?

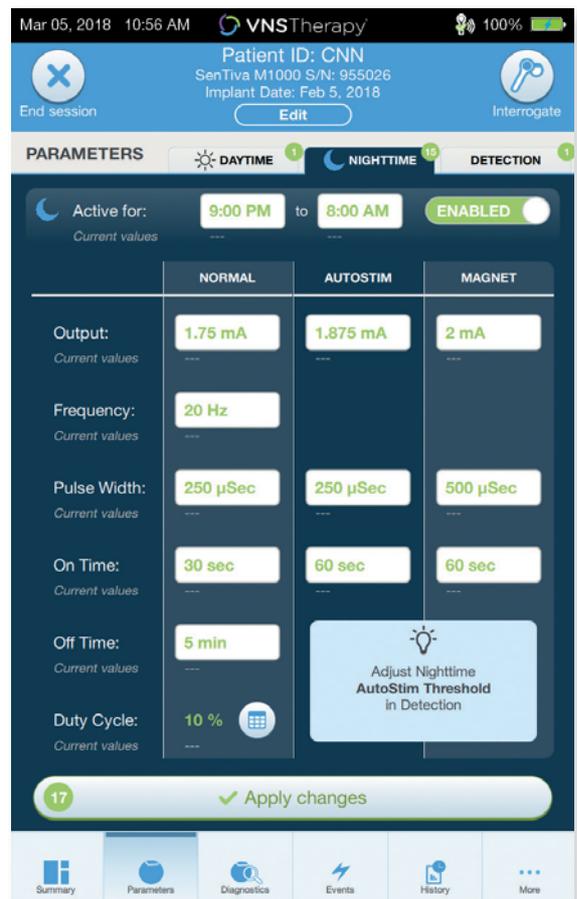

Expanded access to 1.5T and 3T MRI with the latest VNS Therapy technology*2

- The new MRI labeling allows patients with the latest implanted VNS Therapy Systems to receive fast and high quality scans – also with use of Body Coils*


Patients now have access to:



100%
of brain MRI



of all MRI scans
performed on
people with
epilepsy

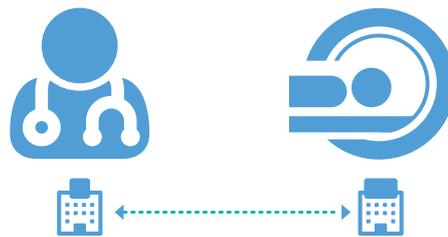

Any MRI center
close by

SenTiva Guidelines for VNS Therapy®

Brief Summary* of Safety Information for the VNS Therapy® System

[Epilepsy Indication] (October 2017)

INTENDED USE / INDICATIONS:

Epilepsy (Non-US)—The VNS Therapy System is indicated for use as an adjunctive therapy in reducing the frequency of seizures in patients whose epileptic disorder is dominated by partial seizures (with or without secondary generalization) or generalized seizures that are refractory to seizure medications. AspireSR® and SenTiva™ feature an Automatic Stimulation Mode which is intended for patients who experience seizures that are associated with cardiac rhythm increases known as ictal tachycardia.

CONTRAINDICATIONS:

Vagotomy—The VNS Therapy System cannot be used in patients after a bilateral or left cervical vagotomy. **Diathermy**—Do not use short-wave diathermy, microwave diathermy, or therapeutic ultrasound diathermy on patients implanted with a VNS Therapy System. Diagnostic ultrasound is not included in this contraindication.

WARNINGS — GENERAL:

Physicians should inform patients about all potential risks and adverse events discussed in the physician's manuals. This document is not intended to serve as a substitute for the complete physician's manuals. The safety and efficacy of the VNS Therapy System have not been established for uses outside the "Intended Use/Indications" chapter of the physician's manuals. The safety and effectiveness of the VNS Therapy System in patients with predisposed dysfunction of cardiac conduction systems (re-entry pathway) have not been established. Post-implant electrocardiograms and Holter monitoring are recommended if clinically indicated. Postoperative bradycardia can occur among patients with certain underlying cardiac arrhythmias. It is important to follow recommended implantation procedures and intraoperative product testing described in the Implantation Procedure chapter of the physician's manual. During the intraoperative System Diagnostics (Lead Test), infrequent incidents of bradycardia and/or asystole have occurred. If asystole, severe bradycardia (heart rate < 40 bpm), or a clinically significant change in heart rate is encountered during a System Diagnostics (Lead Test) or during initiation of stimulation, physicians should be prepared to follow guidelines consistent with Advanced Cardiac Life Support (ACLS). Difficulty swallowing (dysphagia) may occur with active stimulation, and aspiration may result from the increased swallowing difficulties. Patients with pre-existing swallowing difficulties are at greater risk for aspiration. Dyspnea (shortness of breath) may occur with active VNS Therapy. Any patient with underlying pulmonary disease or insufficiency such as chronic obstructive pulmonary disease or asthma may be at increased risk for dyspnea. Patients with obstructive sleep apnea (OSA) may have an increase in apneic events during stimulation. Lowering stimulus frequency or prolonging "OFF" time may prevent exacerbation of OSA. Vagus nerve stimulation may also cause new onset sleep apnea in patients who have not previously been diagnosed with this disorder. Device malfunction could cause painful stimulation or direct current stimulation. Either event could cause nerve damage. Patients should be instructed to use the magnet to stop stimulation if they suspect a malfunction, and then to contact their physician immediately for further evaluation. Patients with the VNS Therapy System or any part of the VNS Therapy System implanted should have MRI procedures performed only as described in the MRI with the VNS Therapy System instructions for use. In some cases, surgery will be required to remove the VNS Therapy System if a scan using a transmit RF body coil is needed. Excessive stimulation at an excess duty cycle (i.e., one that occurs when "ON" time is greater than "OFF" time) and high frequency stimulation (i.e., stimulation at ≥ 50 Hz) has resulted in degenerative nerve damage in laboratory animals. Patients who manipulate the generator and lead through the skin (Twiddler's Syndrome) may damage or disconnect the lead from the generator and/or possibly cause damage to the vagus nerve. Generators with AutoStim only — The AutoStim Mode feature should not be used in patients with clinically meaningful arrhythmias or who are using treatments that interfere with normal intrinsic heart rate responses (e.g., pacemaker dependency, implantable defibrillator, beta adrenergic blocker medications). Patients also should not have a history of chronotropic incompetence (commonly seen in patients with sustained bradycardia (heart rate < 50 bpm)).

WARNINGS — EPILEPSY:

The VNS Therapy System should only be prescribed and monitored by physicians who have specific training and expertise in the management of seizures and the use of this device. It should only be implanted by physicians who are trained in surgery of the carotid sheath and have received specific training in the implantation of this device. The VNS Therapy System is not curative. Physicians should warn patients that the VNS Therapy System is not a cure for epilepsy and that since seizures may occur unexpectedly, patients should consult with a physician before engaging in unsupervised activities, such as driving, swimming, and in strenuous sports that could harm them or others. Sudden unexpected death in epilepsy (SUDEP): Through August 1996, 10 sudden and unexpected deaths (definite, probable, and possible) were recorded among the 1,000 patients implanted and treated with the VNS Therapy device. During this period, these patients had accumulated 2,017 patient-years of exposure. Some of these deaths could represent seizure-related deaths in which the seizure was not observed, at night, for example. This number represents an incidence of 5.0 definite, probable, and possible SUDEP deaths per 1,000 patient-years. Although this rate exceeds that expected in a healthy (nonepileptic) population matched for age and sex, it is within the range of estimates for epilepsy patients not receiving vagus nerve stimulation, ranging from 1.3 SUDEP deaths for the general population of patients with epilepsy, to 3.5 (for definite and probable) for a recently studied antiepileptic drug (AED) clinical trial population similar to the VNS Therapy System clinical cohort, to 9.3 for patients with medically intractable epilepsy who were epilepsy surgery candidates.

PRECAUTIONS — GENERAL:

Physicians should inform patients about all potential risks and adverse events discussed in the VNS Therapy physician's manuals. Prescribing physicians should be experienced in the diagnosis and treatment of depression or epilepsy and should be familiar with the programming and use of the VNS Therapy System.

Physicians who implant the VNS Therapy System should be experienced performing surgery in the carotid sheath and should be trained in the surgical technique relating to implantation of the VNS Therapy System. The safety and effectiveness of the VNS Therapy System have not been established for use during pregnancy. VNS should be used during pregnancy only if clearly needed. The VNS Therapy System is indicated for use only in stimulating the left vagus nerve in the neck area inside the carotid sheath. The VNS Therapy System is indicated for use only in stimulating the left vagus nerve below where the superior and inferior cervical cardiac branches separate from

the vagus nerve. It is important to follow infection control procedures. Infections related to any implanted device are difficult to treat and may require that the device be explanted. The patient should be given antibiotics preoperatively. The surgeon should ensure that all instruments are sterile prior to the procedure. The VNS Therapy System may affect the operation of other implanted devices, such as cardiac pacemakers and implanted defibrillators. Possible effects include sensing problems and inappropriate device responses. If the patient requires concurrent implantable pacemaker, defibrillatory therapy or other types of stimulators, careful programming of each system may be necessary to optimize the patient's benefit from each device. Reversal of lead polarity has been associated with an increased chance of bradycardia in animal studies. It is important that the electrodes are attached to the left vagus nerve in the correct orientation. It is also important to make sure that leads with dual connector pins are correctly inserted (white marker band to + connection) into the generator's lead receptacles. The patient can use a neck brace for the first week to help ensure proper lead stabilization. Do not program the VNS Therapy System to an "ON" or periodic stimulation treatment for at least 14 days after the initial or replacement implantation. For Models 100, 101, 102 and 102R do not use frequencies of 5 Hz or below for long-term stimulation. Resetting the pulse generator turns the device OFF (output current = 0 mA). For Model 100, 101, 102 and 102R resetting the pulse generator will result in device history loss. Patients who smoke may have an increased risk of laryngeal irritation. Generators with AutoStim only — For devices that sense changes in heart rate, false positive detection may cause unintended stimulation. Examples of instances where heart rate may increase include exercise, physical activity, and normal autonomic changes in heart rate, both awake and asleep, etc. Generators with AutoStim only — For the AutoStim feature, the physical location of the device critically affects its ability to properly sense heart beats. Therefore, care must be taken to follow the implant location selection process outlined in the Implantation Procedure. Note that this implant location selection procedure may be performed preoperatively as part of the patient's surgical work-up. M1000 only — Since the Scheduled Programming feature allows the generator to apply therapy increases at scheduled intervals, it may not be appropriate for use in patients who are nonverbal or are unable to use the patient magnet to stop undesired stimulation. Similarly, exercise caution for use of this feature in patients with a history of obstructive sleep apnea, shortness of breath, coughing, swallowing difficulties, or aspiration.

ENVIRONMENTAL AND MEDICAL THERAPY HAZARDS:

Patients should exercise reasonable caution in avoiding devices that generate a strong electric or magnetic field. If a generator ceases operation while in the presence of electromagnetic interference (EMI), moving away from the source may allow it to return to its normal mode of operation. VNS Therapy System operation should always be checked by performing device diagnostics after any of the procedures mentioned in the physician's manuals. For clear imaging, patients may need to be specially positioned for mammography procedures, because of the location of the generator in the chest. Therapeutic radiation may damage the generator's circuitry. Sources of such radiation include therapeutic radiation, cobalt machines, and linear accelerators. The radiation effect is cumulative, with the total dosage determining the extent of damage. The effects of exposure to such radiation can range from a temporary disturbance to permanent damage, and may not be detectable immediately. External defibrillation may damage the generator. Use of electrosurgery [electrocautery or radio frequency (RF) ablation devices] may damage the generator. Magnetic resonance imaging (MRI) should not be performed using a transmit RF body coil for certain VNS Therapy device configurations or under certain specific conditions. In some cases, heating of the lead caused by the transmit RF body coil during MRI may result in serious injury. Static, gradient, and radio frequency (RF) electromagnetic fields associated with MRI may change the generator settings (i.e., reset parameters) or activate the VNS device if the Magnet Mode output remains "ON". Note that certain magnetic resonance (MR) system head coils operate in receive-only mode and require use of the transmit RF body coil. Other MR systems use a transmit/receive RF head coil. Local or surface coils may also be receive-only RF coils that require the transmit RF body coil for MRI. The use of a receive RF coil does not alter hazards of the transmit RF body coil. Exposure of the VNS Therapy System to any transmit RF coil must be avoided. Do not perform MRI scans using any transmit RF coil in the defined exclusion zones. See the MRI with the VNS Therapy System instructions for use for details or further instructions for special cases such as lead breaks or partially explanted VNS Therapy systems. Extracorporeal shockwave lithotripsy may damage the generator. If therapeutic ultrasound therapy is required, avoid positioning the area of the body where the generator is implanted in the water bath or in any other position that would expose it to ultrasound therapy. If that positioning cannot be avoided, program the generator output to 0 mA for the treatment, and then after therapy, reprogram the generator to the original parameters. If the patient receives medical treatment for which electric current is passed through the body (such as from a TENS unit), either the generator should be set to 0 mA or function of the generator should be monitored during initial stages of treatment. Routine therapeutic ultrasound could damage the generator and may be inadvertently concentrated by the device, causing harm to the patient. For complete information related to home occupational environments, cellular phones, other environmental hazards, other devices, and ECG monitors, refer to the physician's manuals.

ADVERSE EVENTS — EPILEPSY

Adverse events reported during clinical studies as statistically significant are listed below in alphabetical order: ataxia (loss of the ability to coordinate muscular movement); dyspepsia (indigestion); dyspnea (difficulty breathing, shortness of breath); hypesthesia (impaired sense of touch); increased coughing; infection; insomnia (inability to sleep); laryngismus (throat, larynx spasms); nausea; pain; paresthesia (prickling of the skin); pharyngitis (inflammation of the pharynx, throat); voice alteration (hoarseness); vomiting. Adverse events reported in clinical investigation of the AutoStim feature were comparable.

*The information contained in this Brief Summary for Physicians represents partial excerpts of important prescribing information taken from the physician's manuals. (Copies of VNS Therapy physician's and patient's manuals are posted at www.livanova.com) The information is not intended to serve as a substitute for a complete and thorough understanding of the material presented in all of the physician's manuals for the VNS Therapy System and its component parts nor does this information represent full disclosure of all pertinent information concerning the use of this product, potential safety complications, or efficacy outcomes.

26-0009-0100/4 (OUS) — 1

SENTIVA™

Fewer seizures.⁵
Shorter seizures.⁶
Faster recovery.^{3,6}
Why wait?



Individual results may vary.
Common side effects include: hoarseness or changes in voice tone,
prickling feeling in the skin, shortness of breath, sore throat and coughing.

LIVANOVA BELGIUM NV
Ikaroslaan 83
1930 Zaventem
Belgium
Tel: +32.2.720.95.93
Fax: +32.2.720.60.53
www.VNSTherapy.com

References:

1. Elliott RE et al. Epilepsy & Behavior 20 (2011) 57-63
2. VNS Therapy Physician's Manual, October 2017, LivaNova, Houston, TX
3. Data on File, LivaNova, Houston, TX
4. Adapted from El Tahry R et al. Epileptic Disord. 2016;18:1-8.
5. Kawai et al, Epileptic Disord 2017; 19 (3): 1-12
6. Orosz et al. Epilepsia 2014, 55(10):1576-1584