

VIOLA™

CLAMPLESS PROXIMAL ANASTOMOSIS SEALING SYSTEM



15.0mm & 17.0mm

CE 0297



Product Description

VIOLA clampless proximal anastomosis sealing system is indicated to maintain hemostasis and facilitate the completion of multiple proximal anastomoses during a coronary artery bypass grafting procedure, without application of an aortic clamp that can lead to neurological complications.


Product Models

Sealing Element Diameter	Model
15.0 mm - 20mm ≤ Aortic Ø < 25mm	FG074
17.0 mm - Aortic Ø ≥ 25mm	FG084

Integrated Accessories

VIOLA PUNCH

A concentric 4.0 mm aortic punch with a detachable handle. The punching element is assembled over the VIOLA catheter wire.

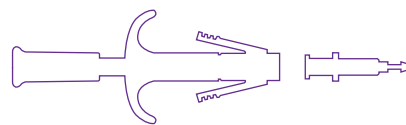

BOUNDARY MARKER

A silicone boundary marker stencil for marking the maximal suture line around the aortic incision.

