


Safire BLU™ Duo Irrigated Ablation Catheter
Therapy™ Cool Path™ Duo Ablation Catheter
Bi-directional Catheter Deflection



ST. JUDE MEDICAL™
MORE CONTROL. LESS RISK.

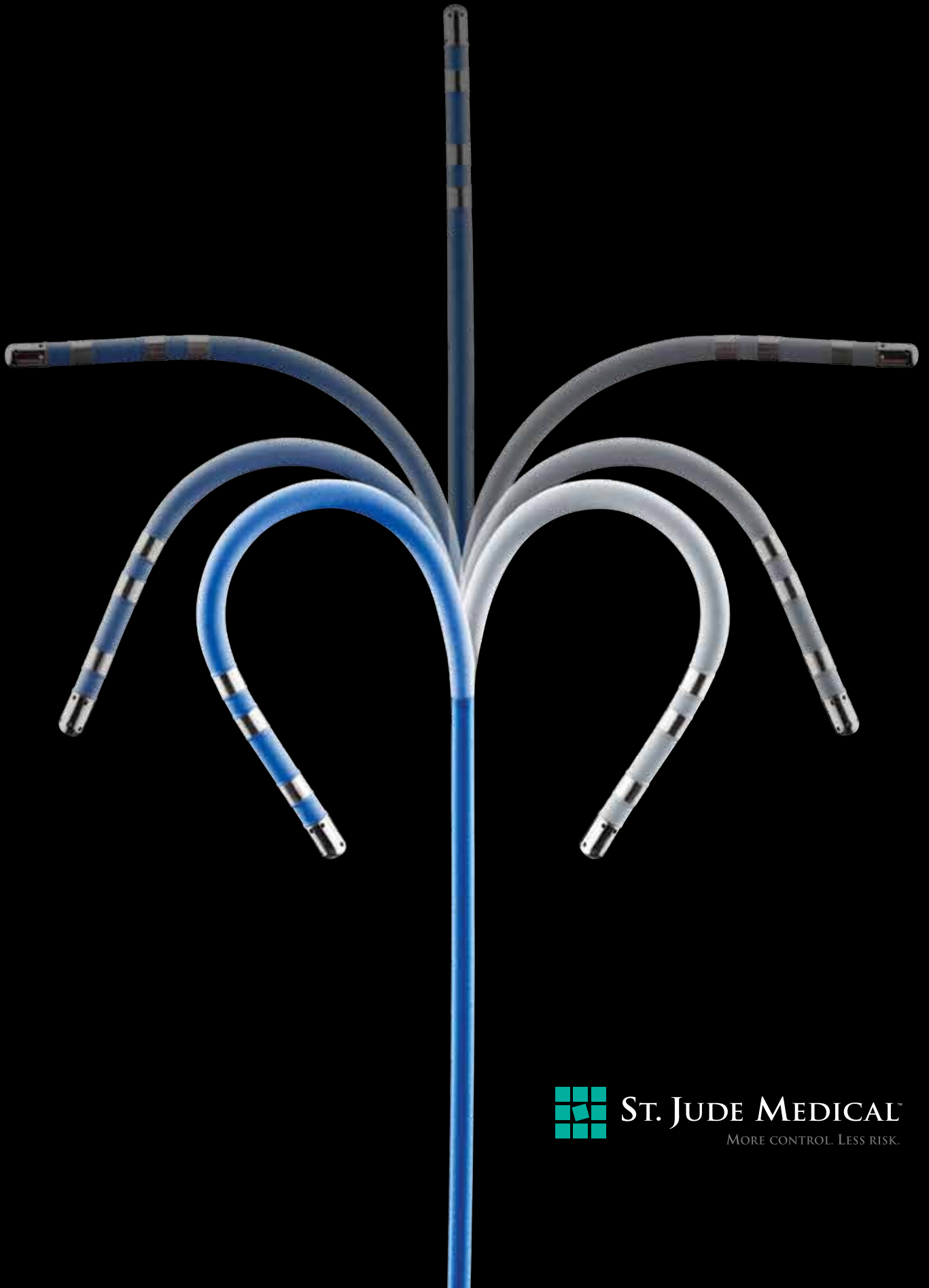

DUO TIP TECHNOLOGY FOR MORE UNIFORM COOLING¹

When the heat is on, Safire BLU™ Duo and Therapy™ Cool Path™ Duo bi-directional ablation catheters feature distal and proximal irrigation ports, providing even cooling across the entire surface of the tip. The high performance shaft with improved torque response offers exceptional control, while bi-directional deflection helps improve reach and maneuverability for easy navigation in challenging anatomies.

High Performance Shaft with Bi-directional Deflection


- New high performance shaft provides improved torque response² offering exceptional catheter control
- Bi-directional catheter deflection provides improved maneuverability to reach challenging anatomic locations during complex procedures

Choice of Handle Tailored for Physician Preference

- Allows physicians to select the style which best suits their preference
 - Safire BLU™ Duo - ComfortGrip™
 - Therapy™ Cool Path™ Duo - Push/pull



Torque Transmission at 180°

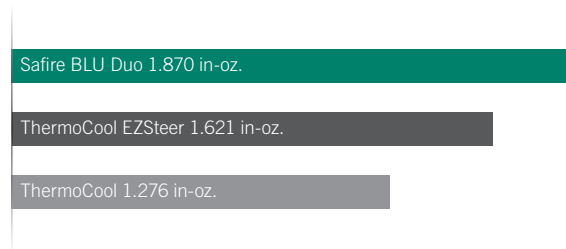

Catheters included in bench testing:

- St. Jude Medical Safire BLU™ Duo irrigated ablation catheter
- Biosense Webster EZSteer™ Thermocool® Bi-directional catheter
- Biosense Webster Thermocool® catheter


- Safire BLU™ Duo catheters had the highest measured torque transmission (1.87 in-oz) at 180° amongst a set of open irrigated catheters during bench testing²
- The Therapy™ Cool Path™ Duo bi-directional ablation catheter uses the same catheter shaft as that of the Safire BLU™ Duo irrigated ablation catheter

Additional Irrigation Ports for More Uniform Temperature Distribution


- Supplemental irrigation ports have been shown to lower the temperature of the proximal area of the electrode¹

- Uniform cooling decreases overall tip temperature¹

Duo Tip Irrigated Ablation Catheter

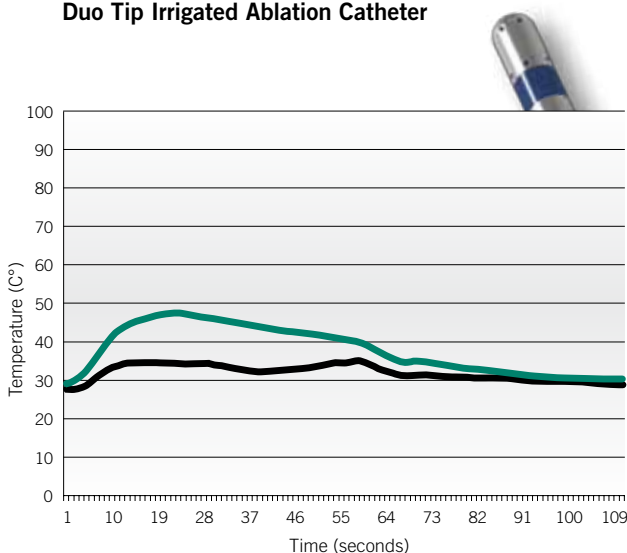

Standard Irrigated Ablation Catheter

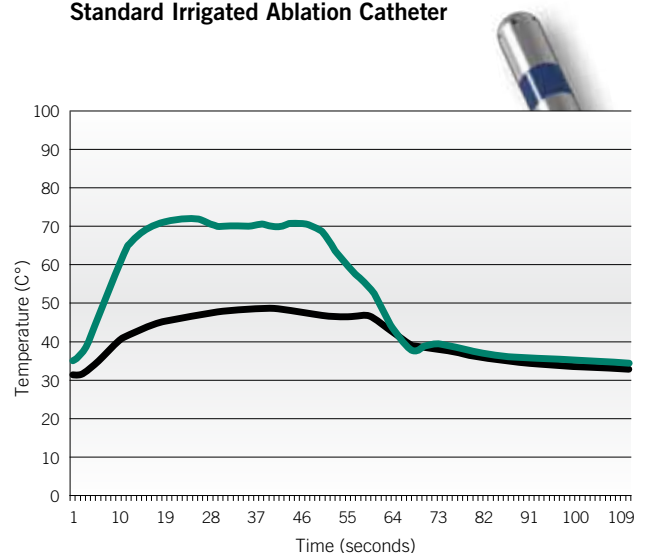

Key Proximal area — Distal area —

1. Experimental conditions: Power = 30 W; flow = 17 ml/min

Ordering Information

Reorder number	Type	Curve	Description	French Size	Tip Size	Electrode Spacing (2 mm rings)	Length
Safire BLU™ Duo Catheter Specifications (Bi-directional)							
A402858	Catheter	Medium Sweep	1304-CP2S-7-25-MS-BD	7 F	4 mm	2-5-2	110 cm
A402859	Catheter	Large Sweep	1304-CP2S-7-25-LS-BD	7 F	4 mm	2-5-2	110 cm
A402860	Catheter	Medium Curl	1304-CP2S-7-25-MC-BD	7 F	4 mm	2-5-2	110 cm
A402861	Catheter	Large Curl	1304-CP2S-7-25-LC-BD	7 F	4 mm	2-5-2	110 cm
Therapy™ Cool Path™ Duo Catheter Specifications (Bi-directional)							
A088050	Catheter	Medium Sweep	1304-CP2-7-25-MS-BD	7 F	4 mm	2-5-2	110 cm
A088052	Catheter	Large Sweep	1304-CP2-7-25-LS-BD	7 F	4 mm	2-5-2	110 cm
A088051	Catheter	Medium Curl	1304-CP2-7-25-MC-BD	7 F	4 mm	2-5-2	110 cm
A088053	Catheter	Large Curl	1304-CP2-7-25-LC-BD	7 F	4 mm	2-5-2	110 cm
Required Catheter Connecting Cables							
85641	Cable	IBI (1641)					250 cm
85713	Cable	Biosense Webster (1713-W)					250 cm

Indications for Use

The Therapy Cool Path Duo and Safire BLU Duo Irrigated Ablation Catheters are intended for creating focal endocardial lesions during cardiac ablation procedures to treat arrhythmias and for cardiac electrophysiology mapping and delivering diagnostic pacing stimuli.

ATRIAL FIBRILLATION

CARDIAC RHYTHM MANAGEMENT

CARDIOVASCULAR

NEUROMODULATION

Global Headquarters

One St. Jude Medical Drive
St. Paul, Minnesota 55117
USA
+1 651 756 2000
+1 651 756 3301 Fax

Atrial Fibrillation Division

One St. Jude Medical Drive
St. Paul, Minnesota 55117
USA
+1 651 756 2000
+1 651 756 3301 Fax

St. Jude Medical Coordination Center BVBA

The Corporate Village
Da Vincilaan 11 Box F1
1935 Zaventem, Belgium
+32 2 774 68 11
+32 2 772 83 84 Fax

St. Jude Medical Brasil, Ltda.

Rua Frei Caneca, 1380
7º ao 9º andares
01307-002 – São Paulo (SP)
Brazil
+55 11 5080 5400
+55 11 5080 5423 Fax

St. Jude Medical (Hong Kong) Ltd.

Suite 1608, 16/F Exchange Tower
33 Wang Chiu Road
Kowloon Bay, Kowloon
Hong Kong SAR
+852 2996 7688
+852 2956 0622 Fax

St. Jude Medical Japan Co., Ltd.

Shiodome City Center 15F
1-5-2, Higashi-Shinbashi
Minato-ku, Tokyo 105-7115
Japan
Tel: +81 3 6255 6370
Fax: +81 3 6255 6371

SJMprofessional.com



ST. JUDE MEDICAL™

MORE CONTROL. LESS RISK.

1 Data on file at St. Jude Medical - Report #90057868

2 Data on file at St. Jude Medical - Report #90056256

Not commercially available in the U.S. Products referenced within are CE marked.

Safire BLU Duo, Therapy Cool Path Duo, ST. JUDE MEDICAL, the nine-squares symbol and MORE CONTROL. LESS RISK. are trademarks and service marks of St. Jude Medical, Inc. and its related companies. ©2010 St. Jude Medical. All Rights Reserved.
Item 100041672EN Rev 01

