

ViewFlex™ PLUS ICE Catheter

HIGH DEFINITION IMAGING, EXCEPTIONAL MANEUVERABILITY WITH SINGLE-HANDED CONTROL

The ViewFlex PLUS Intracardiac Echocardiography (ICE) Catheter is designed for increased range and control of viewing angles. When combined with the market leading Agilis™ handle, tip stability is maximized for enhanced imaging.

Stability, Control and Simplicity

- Large curvature radius provides increased control over viewing angles and positioning
- Steering angles to 120 degrees
- Integrates the market leading Agilis handle for precise single-handed control and full maneuverability during complex procedures
- Exceptional torsional control and tip response


Powerful Integrated Diagnostic Ultrasound System

The ViewFlex PLUS ICE Catheter, together with the ViewMate™ Z Ultrasound Console, provides a full featured, high-fidelity imaging system.

- Unlike fluoroscopy, ICE imaging displays crisp, clear images of anatomical structures.

Imaging Modes

- 2D
- M-Mode
- Pulsed Wave
- Tissue Doppler Imaging
- Continuous Wave
- Color Doppler

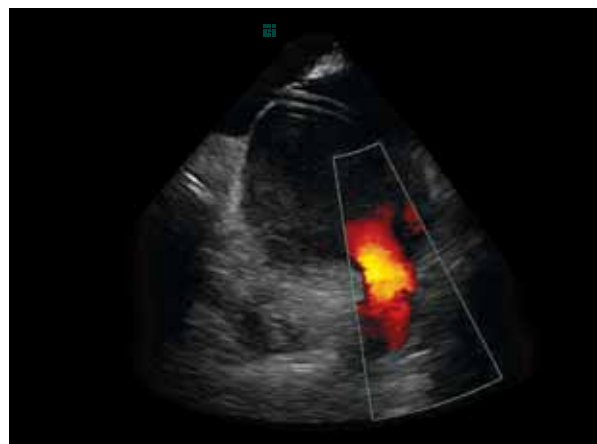

System Safety

- Catheter interface module designed to isolate the patient from the imaging console for patient safety
- Automatic system shut down when catheter temperature threshold is reached
- Integrated sterile connecting cable is designed to prevent potential contaminants from entering the sterile field



ST. JUDE MEDICAL™

MORE CONTROL. LESS RISK.

ViewFlex PLUS ICE Catheter Specifications

Feature	Specification
Catheter size	9 F
Introducer size	10 F
Insertable length	90 cm
Steering	Anterior/Posterior +/- 120 degree flexion
Transducer type	64 element phased array
Viewing sector	90 degrees

Ordering Information

Reorder Number	Product Description
100026714	ViewFlex PLUS ICE Catheter for ViewMate Z Console
100042824	ViewFlex PLUS ICE Catheter for ViewMate II System
100043720	ViewFlex Catheter Interface Module for ViewMate Z Console
100038191	ViewFlex Catheter Interface Module for ViewMate II System

ATRIAL FIBRILLATION

CARDIAC RHYTHM MANAGEMENT

CARDIOVASCULAR

NEUROMODULATION

Global Headquarters

One St. Jude Medical Drive
St. Paul, Minnesota 55117
USA
+1 651 756 2000
+1 651 756 3301 Fax

Atrial Fibrillation Division

One St. Jude Medical Drive
St. Paul, Minnesota 55117
USA
+1 651 756 2000
+1 651 756 3301 Fax

St. Jude Medical

Coordination Center BVBA
The Corporate Village
Da Vincilaan 11 Box F1
1935 Zaventem, Belgium
+32 2 774 68 11
+32 2 772 83 84 Fax

St. Jude Medical Brasil, Ltda.

Rua Frei Caneca, 1380
7° ao 9° andares
01307-002 – São Paulo (SP)
Brazil
+55 11 5080 5400
+55 11 5080 5423 Fax

St. Jude Medical (Hong Kong) Ltd.

Suite 1608, 16/F Exchange Tower
33 Wang Chiu Road
Kowloon Bay, Kowloon
Hong Kong SAR
+852 2996 7688
+852 2956 0622 Fax

St. Jude Medical Japan Co., Ltd.

Shiodome City Center 15F
1-5-2, Higashi-Shinbashi
Minato-ku, Tokyo 105-7115
Japan
Tel: +81 3 6255 6370
Fax: +81 3 6255 6371

SJMprofessional.com



ST. JUDE MEDICAL™

MORE CONTROL. LESS RISK.

Rx Only

Brief Summary: Please review the Instructions for Use prior to using these devices for a complete listing of indications, contraindications, warnings, precautions, potential adverse events and directions for use.

ViewMate, ViewFlex, ST. JUDE MEDICAL, the nine-squares symbol and MORE CONTROL. LESS RISK. are registered and unregistered trademarks and service marks of St. Jude Medical, Inc. and its related companies. ©2011 St. Jude Medical. All Rights Reserved.

Item 100046965EN Rev B

