## ACROSS CTO

## Designed to cross and dilate the toughest lesions.

The **ACROSS CTO<sup>™</sup>** is a PTCA balloon catheter specifically designed for Chronical Total Occlusion. It supports skilled physicians' efforts to provide patients with the benefits of CTO revascularization.

To reach the highest level of performance, the **ACROSS CTO™** PTCA balloon catheter incorporates:

- 1. A unique composite shaft design together with Hydrolubric<sup>™</sup> coating. This results in enhanced pushability and trackability performance in highly tortuous anatomies.
- 2. A truely tapered tip of 4.0 mm in length, to cross safely even the most calcified lesions.

In addition, the **ACROSS CTO™** offers superior high pressure capability due to the balloon material with extreme resistance to potential puncture.

**ACROSS CTO<sup>™</sup>** has currently the **lowest profile on the market**, with a tip profile of 0.015" and a crossing profile of 0.021". It is available with diameters ranging from 1.1 to 2.0 mm, featuring one single marker for more flexibility.


The ACROSS CTO<sup>™</sup> 1.1 is available in Over The Wire (OTW) or in Rapid Exchange (RX) version.



<u>Acrostak</u>

HYDROLUBRIC

## ACROSS CTO™ OTW

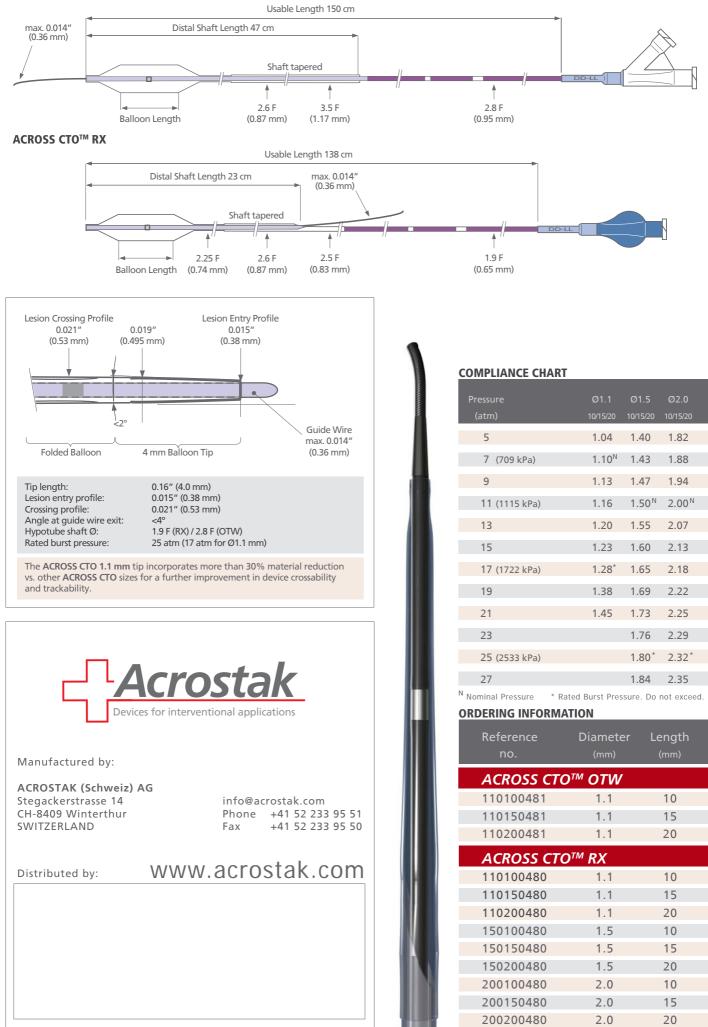

sale in the U.S.

authorised for

Vot