Structural Heart | Valve Repair | Rigid/Semi-Rigid Rings

SJM™ Séguin Semi-Rigid Annuloplasty Ring

Product Highlights


- Designed to aid cardiac output and left ventricular function
- Semi-rigid, three-dimensional posterior section is designed to preserve the physiologic movement of the valve annulus
- Increased rigidity of the anterior section of the ring promotes annular remodeling
- The solid, one-piece inner core resists needle penetration, reducing the potential for suturing through the core

Ordering Information

Contents: Semi-Rigid Ring (1 unit per box)

Reorder Number	Ring Size	Commissure Dimension (mm)	Inside Dimension (mm)	Outside Dimension (mm)	Internal 2-D Orifice Area (mm2)
SARP-24	24	24	22	29	284
SARP-26	26	26	24	31	334
SARP-28	28	28	27	34	404
SARP-30	30	30	28	35	463
SARP-32	32	32	30	37	541
SARP-34	34	34	34	41	602
SARP-36	36	36	36	43	630
SARP-38	38	38	37	45	735
SARP-40	40	40	39	46	815



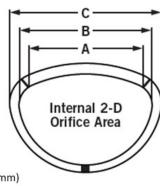

A: Commissure Dimension (mm)

B: Inside Dimension (mm)

C: Outside Dimension (mm)

1. Seguin J, Demaria R, Chaptal PA. Preservation of three-dimensional annular movement with the SJM-Seguin mitral annuloplasty ring. J Heart Valve Dis. 1996;5(6):641-46.

Rx Only

Please review the Instructions for Use prior to using these devices for a complete listing of indications, contraindications, warnings, precautions, potential adverse events and directions for use.

Unless otherwise noted, ® or ™ indicates a registered or unregistered trademark or service mark owned by, or licensed to, St. Jude Medical, Inc. or one of its subsidiaries. ST. JUDE MEDICAL, the nine-squares symbol and MORE CONTROL. LESS RISK. are registered and unregistered trademarks and service marks of St. Jude Medical, Inc. and its related companies. ©2013 St. Jude Medical, Inc. All rights reserved.

100033292 – May 2012

Customer Support: 855-478-5833

Last Updated: April 19, 2013

