

Always next to you

Our team in the UK is always pleased to give you the support needed to optimise your tubing sets. No matter whether you would like to visit our manufacturing facility or make arrangements for one of our dedicated support engineers to visit you. In order to tailor your tubing set and to improve our products and services for you, we do what is most convenient for you. Just contact your local sales representative. Terumo's motto: Always next to you.

Contact Information

Terumo UK • 3 Unity Grove • Off School Lane • Industrial Park South • Knowsley • Merseyside L 34 9 GT • United Kingdom
Tel.: +44 151 545 1573 • Fax: +44 151 545 28 10 • E-mail: coordinatorcx@terumo-europe.com • www.terumo-europe.com/cardiovascular



Terumo Corporation
44-1, 2-chome, Hatagaya
Shibuya-ku, Tokyo 151-0072,
Japan
Tel.: +81 3 3374 8111
Fax: +81 3 3374 8196

European Headquarters

Terumo Europe N.V.
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 39 12 11
Fax: +32 16 40 02 49

European Offices

**Terumo Europe N.V.
Benelux Sales Division**
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 39 25 80
Fax: +32 16 39 25 99
The Netherlands:
Tel.: 0800 0220396
Fax: 0800 0220414

Laboratoires Terumo France S.A.
Parc Ariane - Bât. Uranus
1, rue Hélène Boucher
78284 Guyancourt Cedex, France
Tel.: +33 130 96 13 00
Fax: +33 130 43 60 85

Terumo Deutschland GmbH
Ludwig-Erhard-Strasse 6
65760 Eschborn, Germany
Tel.: +49 6196 80 230
Fax: +49 6196 80 23 200

Terumo Deutschland GmbH
Zweigniederlassung Spreitenbach
Bodenackerstrasse 3
8957 Spreitenbach, Switzerland
Tel.: +41 56 419 10 10
Fax: +41 56 419 10 11

Terumo Europe España SL
Avda. Juan Carlos I, Nº 13 - 7ª Planta
Edificio Torre La Garena
28806 Alcalá de Henares (Madrid),
Spain
Tel.: +34 9021 01 298
Fax: +34 9021 01 358

Terumo Italia S.r.l.
Via Simone Martini, 143/145
00142 Roma RM,
Italy
Tel.: +39 0651 96 14 20
Fax: +39 0650 30 407

Terumo Sweden AB
Nya Varvet, Byggnad 90
Sven Källfelts Gata 18
426 71 Västra Frölunda,
Sweden
Tel.: +46 3174 85 880
Fax: +46 3174 85 890

**Terumo Denmark
Filial of Terumo Sweden AB**
Tel.: +45 7020 93 80
Fax: +45 7020 94 80

Terumo UK Ltd
Tamesis
The Causeway
Egham, Surrey
TW20 9AW, United Kingdom
Tel.: +44 1784 476 200
Fax: +44 1784 476 201

**Terumo Europe N.V.
Moscow Representative Office**
BC Northern Tower
13Fr., 10 Testovskaya Street
Moscow 123317, Russia
Tel.: +7 495 988 4740
Fax: +7 495 988 4739

*Addresses subject to change,
please consult our Website.*

© Registered Trademark
Published by Terumo Europe N.V.
CV147GB-0811FKPR-I(0811)E

Cardiovascular Tubing Sets

European Engineering and Tubing Set Production



Adapted to your needs

European Engineering and Tubing Set Production

Terumo is committed to manufacturing quality perfusion tubing sets for cardiac surgery. The Terumo Europe N.V. UK production site, located between Manchester and Liverpool, was opened in 1996 and has since increased its capacity enormously. Today the Terumo UK facility produces over 64,000 tubing sets annually and coordinates the work of external tubing set suppliers throughout Europe.



The cardiovascular healthcare community demands products they can depend on. Terumo delivers quality – dependably, predictably and consistently.

The Terumo UK facility's range of perfusion circuit products is one of the most comprehensive in the industry, offering solutions from the smallest circuit to a fully pre-connected oxygenator tubing set. At the Terumo UK plant, design engineers work directly with customers to create tailored solutions. From drawing to validation and production, Terumo professionals provide technical support focused on helping perfusion teams to give the very best patient care.

Our Products

- **Tailored Tubing Sets**
- **Standard Products:**
 - Haemoconcentrator Sets
 - Cardioplegia Systems
 - Accessory Sets
- **Modular Perfusion Circuits:**
 - ROCSafe
 - Prescriptive Oxygenation



Quality

Terumo's state-of-the-art manufacturing facilities in UK make extensive use of modern processes and systems. Our quality system is certified to ISO 9001 and ISO 13485 standards. Terumo's superior manufacturing processes are a vivid example of the success that can be attained when people focus on quality. We apply the strictest standards of quality and control to every item we manufacture.



Xcoating is Terumo's exclusive biocompatible surface coating. Xcoating is composed of an amphiphilic polymer, which means that it has both hydrophobic and hydrophilic properties. Xcoating is biopassive and will not react with blood components. It is not heparin-based and can be used with heparin-intolerant patients. Xcoating can be applied to virtually all device material surfaces, including polypropylene, polycarbonate, polyurethane, stainless steel, filter mesh, and PVC.



Terumo has established strict guidelines to ensure people and environmentally friendly product development. We promote mercury-free, PVC-free and DEHP-free products and we develop energy saving models for devices that are left on for long periods of time.




TOTM – an alternative plasticizer
Terumo is ever striving to develop new medical technologies with minimal burden to patients and the environment. Our efforts are closely linked to our corporate philosophy of contributing to society through health care by providing safe and trustworthy products. In line with this vision, Terumo provides an alternative plasticizer for the manufacturing of its products. TOTM (trioctyl trimellitate) offers outstanding physical properties, such as flexibility, to the material and low plasticizer elution.

Packaging configurations


Pre-connected Tray Set

Tailored complete tubing set with oxygenator provided in one Tyvek lidded tray. This packing style accommodates a fully pre-connected tubing set reducing set up time which can be crucial in emergency cases.

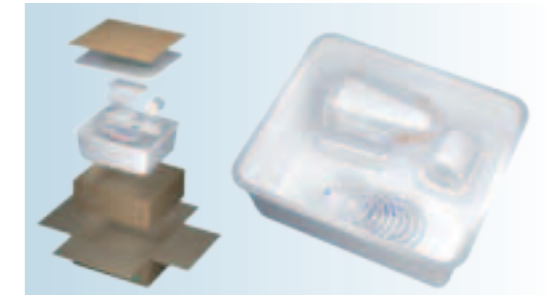


Masterbox

Tubing set and oxygenator packed next to each other in one box. Environmental friendly option requires less packaging. The packaging material is biodegradable.

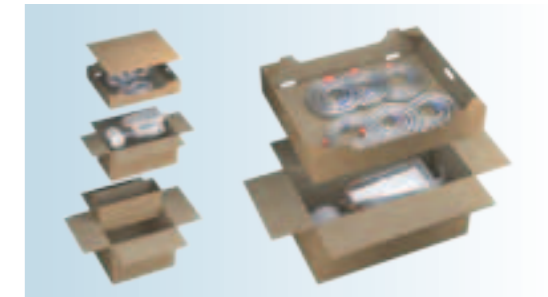


Yellow Bin

Tubing set and oxygenator provided in a disposable bin. Everything you need in one box cuts down on the number of orders needed. Less storage space required. Can be provided partially pre-connected for easier set up. Reduces packaging waste and fulfills all requirements of waste incinerating.

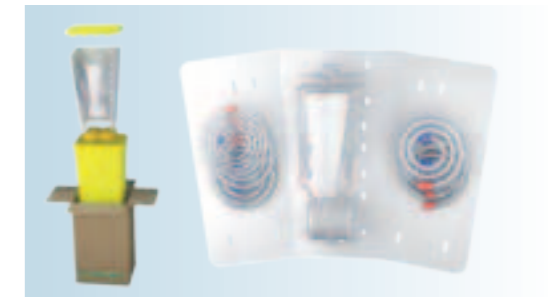

Eco Packaging

Tailored tubing set provided in our most environmentally friendly packaging featuring easy disposal and a smaller carbon footprint. Using minimal packaging without compromising sterility and quality, our sets range from complex to standard.

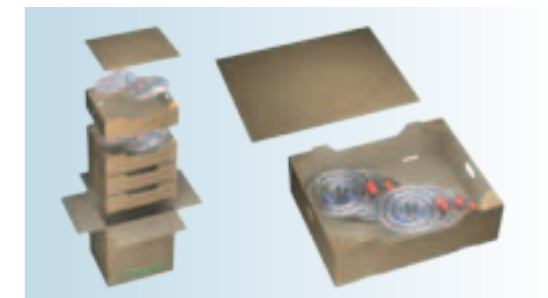

Tubing Set Tray

Tailored tubing set in Tyvek lidded tray. Sets range from complex to standard with a variety of tray sizes available allowing complete flexibility in tubing set design.

